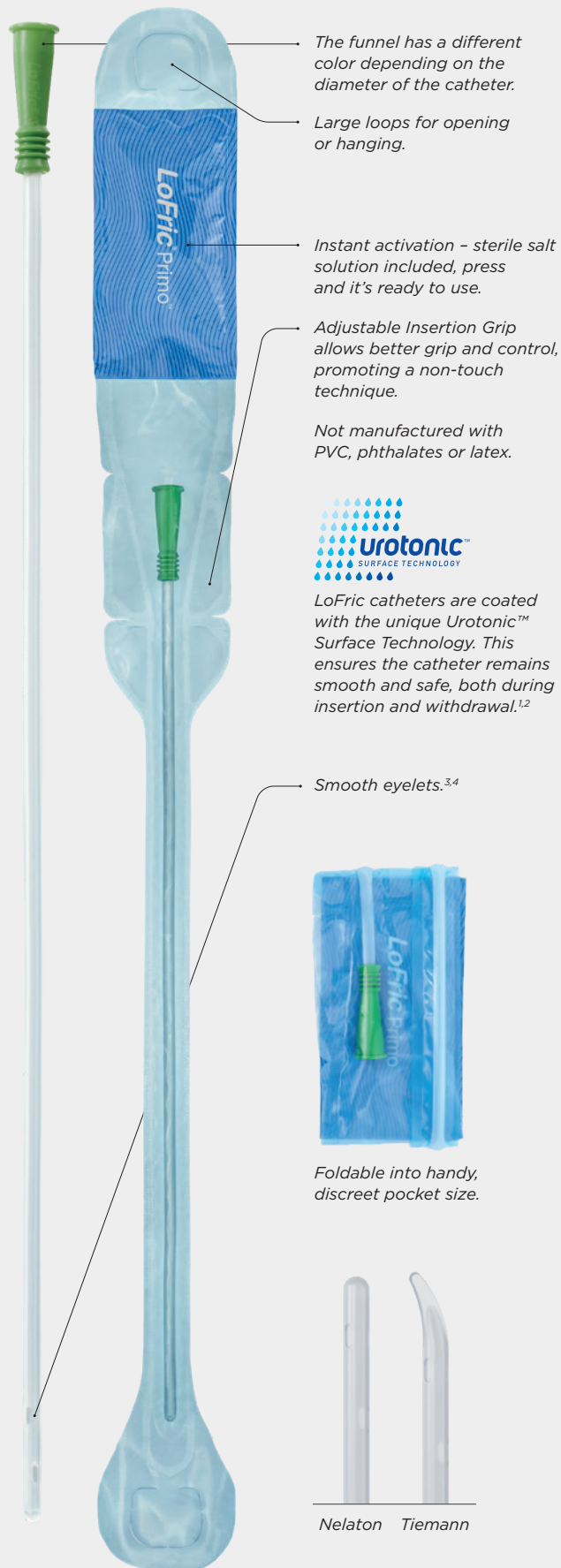




The funnel has a different color depending on the diameter of the catheter.

Large loops for opening or hanging.

Instant activation - sterile salt solution included, press and it's ready to use.

Adjustable Insertion Grip allows better grip and control, promoting a non-touch technique.

Not manufactured with PVC, phthalates or latex.



LoFric catheters are coated with the unique Urotonic™ Surface Technology. This ensures the catheter remains smooth and safe, both during insertion and withdrawal.^{1,2}

Smooth eyelets.^{3,4}



Foldable into handy, discreet pocket size.



Nelaton Tiemann

Product description

LoFric® Primo™, single use, hydrophilic urinary catheter with salt solution for activation, intended for intermittent urinary catheterization.

Sterilization and shelf life

Sterilized by irradiation. Use before expiry date on package.

Package

30 catheters per customer box. 2 or 4 customer boxes per shipping box.

Material composition

Hydrophilic surface coating: Polyvinyl pyrrolidone (PVP).

Catheter: Polyolefin-based elastomer (POBE).

Sterile solution pocket: Laminate foil of polyethylene terephthalate (PET), aluminum and polyethylene (PE). Sterile water with sodium chloride (NaCl).

Item packaging: Polyethylene (PE).

Product range

Reference	Type	Length	Size
410 08 XX*	LoFric® Primo™ Nelaton (male)	40 cm	CH 08
410 10 XX*	LoFric® Primo™ Nelaton (male)	40 cm	CH 10
410 12 XX*	LoFric® Primo™ Nelaton (male)	40 cm	CH 12
410 14 XX*	LoFric® Primo™ Nelaton (male)	40 cm	CH 14
410 16 XX*	LoFric® Primo™ Nelaton (male)	40 cm	CH 16
410 18 XX*	LoFric® Primo™ Nelaton (male)	40 cm	CH 18
411 06 XX*	LoFric® Primo™ Nelaton (paediatric)	20 cm	CH 06
411 08 XX*	LoFric® Primo™ Nelaton (paediatric)	20 cm	CH 08
411 10 XX*	LoFric® Primo™ Nelaton (paediatric)	20 cm	CH 10
413 08 XX*	LoFric® Primo™ Nelaton (female)	20 cm	CH 08
413 10 XX*	LoFric® Primo™ Nelaton (female)	20 cm	CH 10
413 12 XX*	LoFric® Primo™ Nelaton (female)	20 cm	CH 12
413 14 XX*	LoFric® Primo™ Nelaton (female)	20 cm	CH 14
413 16 XX*	LoFric® Primo™ Nelaton (female)	20 cm	CH 16
413 18 XX*	LoFric® Primo™ Nelaton (female)	20 cm	CH 18
414 08 XX*	LoFric® Primo™ Nelaton (female)	15 cm	CH 08
414 10 XX*	LoFric® Primo™ Nelaton (female)	15 cm	CH 10
414 12 XX*	LoFric® Primo™ Nelaton (female)	15 cm	CH 12
414 14 XX*	LoFric® Primo™ Nelaton (female)	15 cm	CH 14
415 10 XX*	LoFric® Primo™ Tiemann (male)	40 cm	CH 10
415 12 XX*	LoFric® Primo™ Tiemann (male)	40 cm	CH 12
415 14 XX*	LoFric® Primo™ Tiemann (male)	40 cm	CH 14
415 16 XX*	LoFric® Primo™ Tiemann (male)	40 cm	CH 16
415 18 XX*	LoFric® Primo™ Tiemann (male)	40 cm	CH 18

*Generic article number without the two digit suffix specific for a region or country destination when distributing an article.



User-friendly packaging.



LoFric® Primo™

References

1. Waller et al., Spinal Cord, 1997
2. Dunger et al., Acta Endocrinologica, 1988
3. Lundgren et al., Spinal Cord, 2000
4. Newman et al., Springer, 2018

At Wellspect we develop innovative continence care solutions that improve quality of life for people with bladder and bowel problems. We inspire our users to build self-confidence and independence as well as good health and well-being. We have been leading the industry for over 40 years with our product brands LoFric[®] and Navina[™]. We always aim to minimize the environmental impact of our products and passionately strive to become climate neutral. We work together with users and healthcare professionals to improve clinical outcome in a sustainable way, now and for the future.

Wellspect. A Real Difference.

For more information about our products and services, please visit [Wellspect.com](https://www.wellspect.com).

Join the conversation on LinkedIn, Twitter, Facebook and Instagram.

[wellspect.com](https://www.wellspect.com)

Wellspect HealthCare, Aminogatan 1, P.O. Box 14, SE-431 21 Mölndal, Sweden.
Phone: +46 31 376 40 00. www.wellspect.com

LoFric



Manufacturer: Wellspect HealthCare, Aminogatan 1, P.O. Box 14, SE-431 21 Mölndal, Sweden.
Phone: +46 31 376 40 00. www.wellspect.com

CE 2797

R_x Only